A E N O



Quick Start Guide

Toothbrush DB7 / DB8


(Models: ADB0007,
ADB0008)

v 2.2.2 | 03.09.2025


ENG.....	6
ARA.....	8
BOS.....	10
BUL.....	12
CES.....	14
DEU.....	16
ELL.....	18
EST.....	20
FRA.....	22
HRV.....	24
HUN.....	26
HYE.....	28
ITA.....	30
KAT.....	32
KAZ.....	34
LAV.....	36
LIT.....	38
NLD.....	40
POL.....	42
POR.....	44
RON.....	46
RUS.....	48
SLK.....	50
SLV.....	52
SPA.....	54
SRP.....	56
UKR.....	58
UZB.....	60

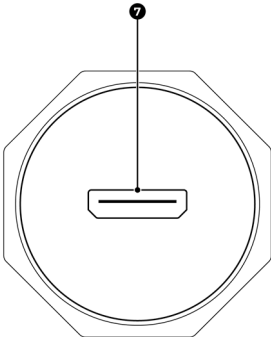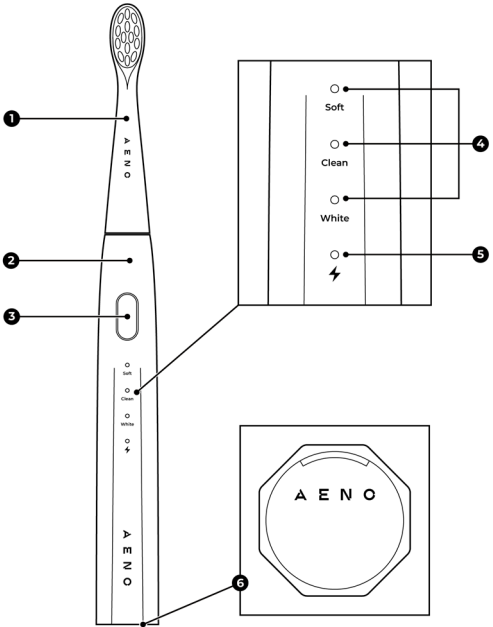
B

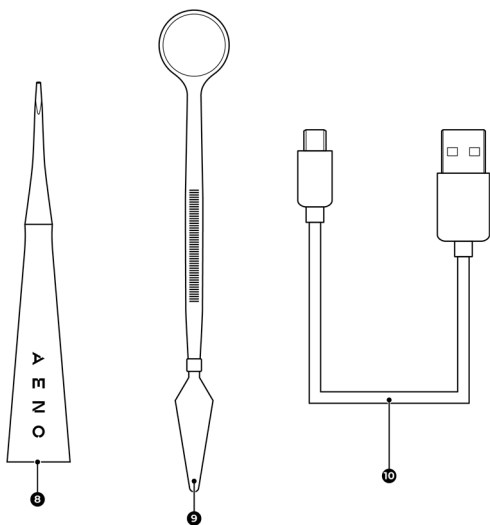
C

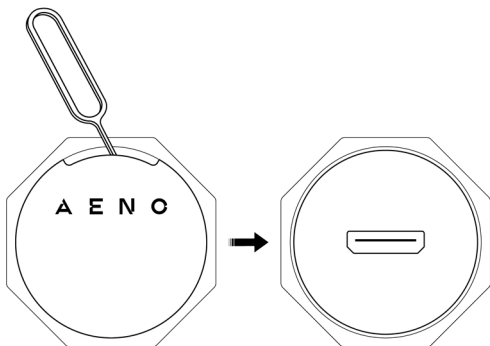D

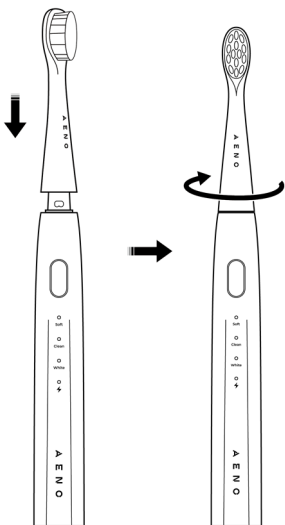E

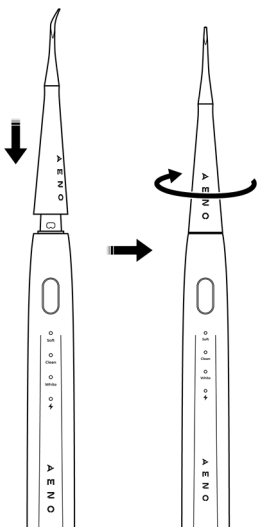ENG AENO™ Toothbrush is intended for daily oral hygiene.

Technical Specifications

Power supply (input): 5.0 V / 0.3 A (DC). Battery: lithium-ion, non-removable, 3.7 V, 800 mAh. Full charge time: up to 3 h. Operating time on a single charge: up to 10 h. Ingress protection rating: IPX7. Number of operating modes: 3. Vibration intensity: up to 30,000 vibrations/minute. Bristle stiffness: medium. Service life of brush heads: 3 months. Noise level (at a distance of 15 cm): ≥ 55 dB. Color: White. Control: manual. Casing material: ABS plastic, TPE. Toothbrush dimensions without the brush head (L×W×H): 21×21×163 mm. Weight (without the brush head): 49 g. USB (A) / micro-USB cable length: 0.3 m. Operating conditions: temp. 0...+40 °C, RH 0–95 %. Storage conditions: temp. 0...+40 °C, RH 0–75 %.

Package Contents (see Fig. A, B)

ADB0007: Toothbrush, replacement brush head (1 pc.), USB (A) / micro-USB cable, sticker sheets (2 pcs.), Quick Start Guide.

ADB0008: Toothbrush, replacement brush heads (3 pcs.), plaque remover, dental mirror, USB (A) / micro-USB cable, Quick Start Guide.

Device Elements and Accessories (see Fig. A, B)

1 – replacement brush head, 2 – handle, 3 – power button, 4 – operating mode indicators, 5 – battery charge indicator, 6 – protective cover, 7 – micro-USB port, 8 – plaque remover, 9 – dental mirror, 10 – USB (A) / micro-USB cable.

Restrictions and Warnings

This device is not intended for use by children under the age of 6 or by persons with reduced physical, mental or intellectual abilities unless they have sufficient experience and knowledge of how to use the device and are supervised by a person responsible for their safety. Do not allow children to play with the device and its accessories. Do not drop or throw the device. Do not turn on or use a damaged device. If you have a pacemaker, oral diseases, or if you have just undergone oral surgery, consult your doctor before using the device. If your gums begin to bleed after using the device, and the bleeding continues for several days, contact your dentist. Each family member should use an individual brush head. Replace the brush head at least once every 3 months. Use only original replacement heads. The plaque remover is intended for use by adults only. If you are not using the device for a long time, charge it every 2–3 months to prevent deep discharge of the battery. Do not wash the device or its accessories in a dishwasher. Do not use chemical and aggressive detergents, abrasive pastes, products containing acids and solvents or metal sponges for cleaning. For a detailed description of the device, its modes and features, see the complete User Manual available at aeno.com/documents.

Statuses of the Battery Charge Indicator

Color	Status	Meaning
Green	On	The device is turned on. Battery level is medium/high
		The device is connected to a power supply. Battery is fully charged
Yellow	On	The device is connected to a power supply. Charging is in progress
Red	On	The device is turned on. Battery level is low
	Blinking	The battery is fully discharged
–	Off	The device is not connected to a power supply or is turned off

Operating Modes Indicators

When you select the operating mode, the corresponding indicator (4) lights up white.

Power Button Functions

Action	Meaning
Press once (device is turned off)	Turn the device on
Press once (the device is turned on, and it has been less than 5 seconds since it was turned on)	Switch operating modes
Press once (the device is turned on, and it has been more than 5 seconds since it was turned on)	Turn the device off

If you have any questions or difficulties in using the AENO™ device, please e-mail support at support@aeno.com or chat with us online at aeno.com/service-and-warranty. The specialists will help you solve the problem, and you will not have to waste time and effort visiting the store.

Preparing for Use

Carefully open the packaging and remove the device and its accessories. Clean them (see "Cleaning and Maintenance"). Remove the protective cover (6) from the brush handle (see Fig. C).

ATTENTION! The protective cover ensures that the micro-USB port (7) is waterproof. If the protective cover is removed, do not allow liquids to enter the port of the device.

Note. A power adapter is not included. The adapter output parameters must be 5.0 V / 1.0 A (DC).

Connect the device to a power source using a USB (A) / micro-USB cable (10). The battery charge indicator (5) will light up yellow. When the indicator turns green, the battery is fully charged. Disconnect the device from the power source, and then disconnect the cable from the micro-USB port. Insert the protective cover into the device handle.

Choosing Operating Modes of the Device

Turn on the device by pressing the power button (3) once. Subsequent presses of the button will switch the operating modes. The last used mode is saved in the device settings the next time the device is turned on.

Operating Modes

- "Soft". Delicate, gentle brushing of the teeth. Recommended for brushing sensitive teeth, as well as for brushing children's teeth.
- "Clean". Daily brushing. Recommended for standard brushing.
- "White". Teeth whitening. Recommended for people who smoke and drink tea or coffee often.

Use

With the brush head. Install the head (1) as shown in Fig. D. Wet the bristles of the brush head and apply the desired amount of a dental care product to the bristles. Point the bristles toward the surface of the teeth at an angle of about 45°. Press the power button to turn on the device. Brush your teeth without putting too much pressure on the enamel. Every 30 seconds the device will notify you to change the brushing zone by briefly stopping the vibration. Change the brushing zone after each notification. After 2 minutes of operation, the device will turn off automatically.

With plaque remover. Install the attachment (8) as shown in Fig. E. Point the metal end of the attachment toward the tooth surface at an angle of about 45°. Press the power button to turn on the device. Clean problem areas using gentle vertical strokes.

ATTENTION! Do not apply excessive force when using the plaque remover. Careless use of the attachment can cause oral trauma as well as damage tooth enamel.

Note. Using a plaque remover cannot substitute professional oral hygiene at the dentist.

Cleaning and Maintenance

ATTENTION! Do not clean the device until it has been turned off and disconnected from the power source.

Rinse the device and accessories under running water. Then wipe dry.

Troubleshooting

The toothbrush does not turn on. Possible cause: the battery is discharged. Solution: charge the device (charging time – at least 3 hours).

The device is connected to a power source but is not charging (the battery charge indicator is not on). Possible cause: the device has not been used for a long time. Solution: connect the toothbrush to a power source and wait about 15 minutes. If the battery charge indicator does not light up, contact an authorized service center.

ATTENTION! If none of the suggested solutions helped you eliminate the problem, contact the supplier or an authorized service center. Do not disassemble the device or attempt to repair it yourself.

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ENG ASBISC reserves the right to modify the device and to revise and change this document without prior notice to users.

The warranty period and service life are 2 years from the date of retail sale of the device. Manufacturer: ASBISC Enterprises PLC, Iapetou, 1, Agios Athanasios, 4101, Limassol, Cyprus. Tel.: +48732080077, info@asbis.com, asbis.com. Made in China. All trademarks used herein are the property of their respective owners. The date of manufacture is marked on the package.

Up-to-date information about and a detailed description of the device, as well as connection instructions, certificates, and details about the companies that deal with quality and warranty claims are available at aeno.com/documents. This device complies with applicable safety requirements.

RoHS Directive. The device complies with the requirements of the Directive 2011/65/EU on the Restriction of Hazardous Substances (RoHS), including the requirements of the RoHS Directive 2015/863/EU.

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ENG Recycling information. These symbols indicate that you must follow the Waste Electrical and Electronic Equipment (WEEE) and Waste Battery and Battery Regulations when disposing of the device, its batteries and accumulators, and its electrical and electronic accessories. According to the rules, this equipment must be disposed of separately at the end of its service life. Do not dispose of the device, its batteries and accumulators, or its electrical and electronic accessories together with unsorted municipal waste, as this will harm the environment. To dispose of this equipment, it must be returned to the point of sale or turned in to a local recycling facility. You should contact your local household waste disposal service for details.

For devices with a non-removable battery: under normal or abnormal operating conditions, as well as in the event of a single malfunction, the charging voltage and current must not exceed their maximum permissible values. If the battery temperature exceeds the maximum allowable temperature, then, under abnormal conditions, the device will stop charging. If the battery temperature drops below the minimum allowable temperature during charging, then, under abnormal conditions, the current will be limited.

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